


Duette[®]

Multi-Band Mucosectomy

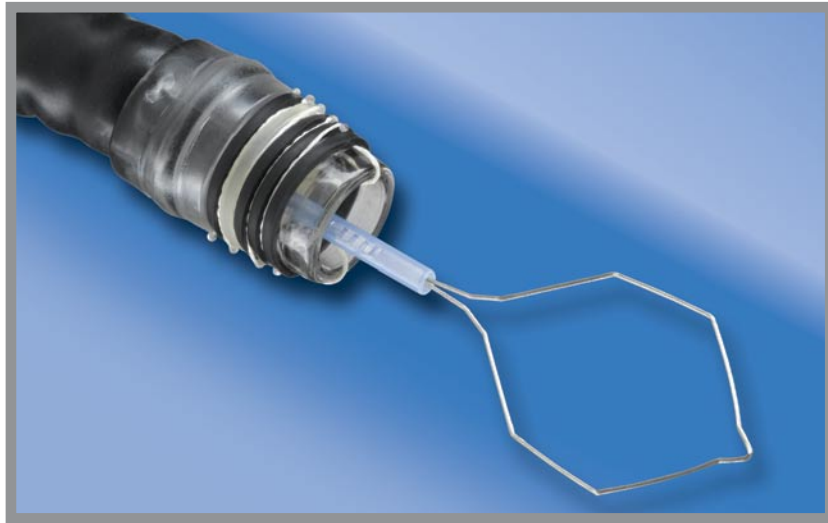

Cook Endoscopy presents Duette, an innovative solution for endoscopic mucosal resection. Duette is a combined multi-band ligator device and snare which allows simple ligation and snare resection of superficial lesions and early cancers in the upper GI tract.

COOK[®]
Endoscopy



DUETTE®

MULTI-BAND MUCOSECTOMY

Duette is intended for endoscopic mucosal resection in the upper GI tract. This device is intended for single use only. With its integrated components Duette offers efficiency and simplicity during endoscopic mucosal resection. The snare component is supplied sterile and the bander component is supplied non-sterile.

- ▶ Easy set-up, similar to Saeed Six-Shooter® Multi-Band Ligator
- ▶ Six latex rubber bands for resection of large lesions
- ▶ Soft, mini-hexagonal, braided snare cuts cleanly and has excellent shape retention even with repeated resections
- ▶ Snare may be passed through ligator handle allowing for speedy resection of ligated tissue without removal of endoscope or disassembly of components
- ▶ Minimum endoscope channel diameter 3.7 mm



Order Number	Ligator		Snare		Minimum Accessory Channel
	Endoscope Outer Diameter	Number of Bands	Sheath Size	Size	
DT-6	9.5 - 13 mm	6	7.0 FR	1.5 x 2.5 cm	3.7 mm
DT-6-XL	11 - 14 mm	6	7.0 FR	1.5 x 2.5 cm	3.7 mm

COOK®
Endoscopy

US Distributor: Cook Endoscopy
4900 Bethania Station Road, Winston-Salem, NC 27105
www.cookendoscopy.com

Duette and Saeed Six-Shooter are registered trademarks of Wilson-Cook Medical Inc.
Cook is a registered trademark of Cook Incorporated.